


NORLAKE[®] MINI-ROOM WALK-INS

TEN DAY SHIPMENT

Controlled Temperature Chambers-Stability Rooms

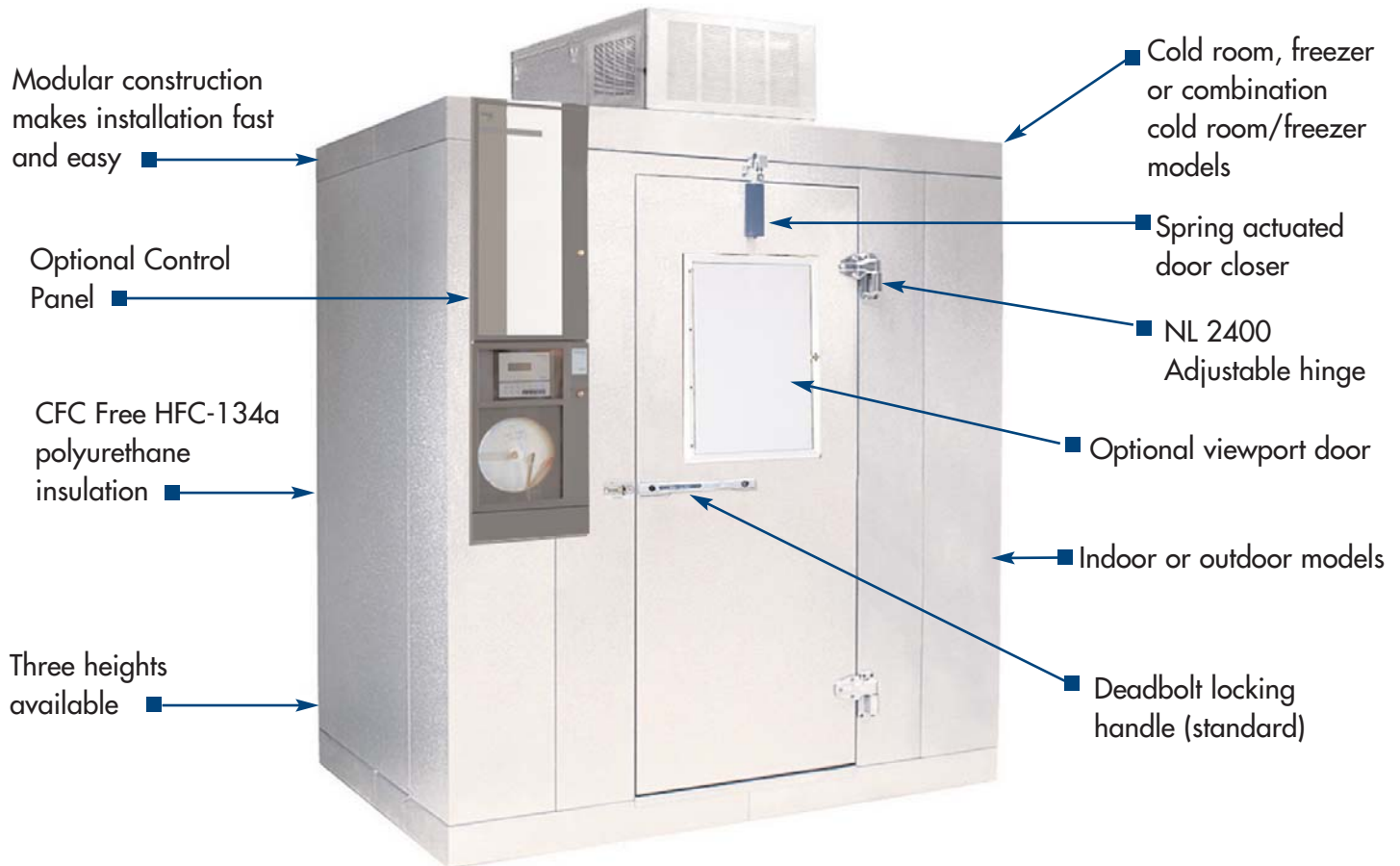

10,000 Models Available to Ship in Ten Days!

Standard Features

- Sizes from 4' x 5' to 8' x 14'
- Three Heights: 6'7" & 7'7" with floor, 7'4" floorless
- Temperatures: +4°C, -20°C, -30°C, 4°C to 55°C
- Exterior finish: Corrosion resistant stucco embossed temper passed coated steel
- 26" wide self closing door
- 15 year panel warranty

Options Available for 10 Day Shipment

- Microprocessor based control panel with accuracies to +0.15°C
- Humidification
- Dehumidification
- Chart recorder
- Viewport door
- Duplex electrical outlets (surface mounted)
- Hi/Lo alarm
- Exterior ramp
- Shelving kits
- Heater kits

Mini-Room Warm Rooms

- Mini-Room Warm Rooms (4°C to 55°C)
- Mini-Room Warm Rooms are equipped with a Nor-Lake Scientific Control Panel and a ceiling mount refrigeration system.

Call 800-477-5253 Today for a Mini-Room or a Custom Enviro Line™ Room Designed Specifically for You.

NORLAKE[®] SCIENTIFIC

NOR-LAKE, SCIENTIFIC
727 Second Street
Hudson, Wisconsin 54016
800-477-5253
715-386-6149 FAX
www.norlake.com

MINI-ROOM WALK-INS • CONTROLLED TEMPERATURE CHAMBERS

Multiple Applications

- Pharmaceutical Research and Development
- Drug Stability and Shelf Life Testing
- Plasma and Blood Storage
- Clinical Research
- Environment Simulation Testing
- Working Laboratories
- Biological Research
- Electronic Stress Testing

Mini-Room Walk-in available as...

- Cold room, freezer or combination cold room/freezer models
- Indoor or outdoor models
- Sizes from 4' x 5' to 8' x 14'
- Three Heights:
6'7" & 7'7" with floor, 7'4" floorless
- Capacity from 77 cubic feet to 1,000 cubic feet
- Wall or ceiling mounted Capsule Pak™ Refrigeration Systems
- Temperatures: +4°C, -20°C, -30°C, 4°C to 55°C

The Mini-Room Walk-in offers you...

- Installation that is fast and easy
- A full 4" panel foamed-in-place with CFC and HCFC Free HFC-134a polyurethane insulation
- Deadbolt locking handle for security
- Nor-Lake Scientific's exclusive Capsule Pak™ refrigeration system that features a flush coil. 100% usable interior space
- Automatic condensate vaporizer. No plumbing required on indoor units
- New Sizes: Expanded Mini-Room Cold Rooms
 - 6' x 14'
 - 8' x 12'
 - 8' x 14'
 - 8' x 20' combination
- Smaller Mini-Room Walk-ins
 - 3'6" x 6'0" x 6'7"
 - and
 - 3'6" x 7'0" x 6'7"

Similar floor space as 3 door laboratory reach-in with 20%-30% more capacity

Capsule Pak™ Refrigeration Systems



Capsule Pak™ Refrigeration Systems consist of a unitized system which is factory assembled, wired, charged, tested and fully equipped for insertion into a factory prepared walk-in opening. The Capsule Pak™ System has a flush coil which keeps all components outside, allowing more storage inside.


- Wall or ceiling mounted systems
- Indoor or outdoor models
- Temperatures: 4°C, -20°C, -30°C, 4°C to 55°C
- Optional Remote Capsule Pak™ Systems
- U.L. & C.S.A. listed



Shelving Kits

Specially designed Shelving Kits are available for all Mini-Room walk-ins.

- Quick, easy assembly
- Rigid construction
- Stronger shelves
- Flexible shelf adjustment-2" increments
- By using the Nor-Lake Klip™ the shelving layouts utilize maximum space for storage leaving the door opening clear of end standards
- Green Kote™ special epoxy coated finish provides extra protection from acidic materials

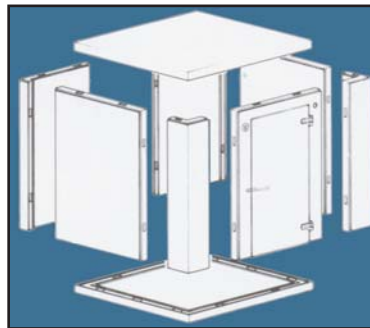

Modular Assembly



ISO 9001:2000
Registered to ISO 9001:2000
File No. A3204

Standex
Food Service

NORLAKE SCIENTIFIC

727 Second Street
Hudson, Wisconsin 54016
800-477-5253
715-386-2323
Revision Date: 11/05
© 2005 NOR-LAKE, INC.
Part Number: 089730