

# STAINLESS STEEL DELUXE STACKING CAT CAGES WITH REMOVABLE COLLECTION PANS

Series A Cages fit Series A Racks



Standard Full Depth Floor



Optional 1/2 Floor and Litter Pan

## 101447-00 Cage (shown)

This Stainless Steel Deluxe Stacking Cat Cage with Removable Collection Pan features a one-piece wrapper, consisting of one solid side, solid back, and solid lower portion of the opposite side.

The top and upper portion of one side of the cage have electropolished grills. The wrapper and grills are welded to the cage's welded angle frame.

All welds are polished to a satin finish, removing any discoloration.

Each cage is equipped with either one or two shelves, depending on the cage's height. Shelves are available in three types of materials: stainless steel, oak, or polyethylene.

The removable electropolished door has a rod-frame with vertical and horizontal rods, a 5" wide by 3" high cardholder, a rotating lazy-susan feeder with two removable 12 ounce bowls, and a two-point self-locking latch.

The removable electropolished floor with rolled-edge strip animal resting surface, has a frame design that holds the floor off of the wrapper's deflector. This allows waste to drain

unobstructed into the removable tapered collection pan. The collection pan is supported by the cage's pan runners below the deflector drip edge. This unique drip edge directs waste into the collection pan.

The cage's lower front crossmember is designed to hold the cage securely on the rack.

### OPTIONS

603672 Stainless Steel Litter Pan 13"L 9"W 4"D

#### Stainless steel rolled edge 1/2 floors\*


207852 For 3 sq ft cage  
207853 For 4 sq ft cage  
207854 For 5 sq ft cage

#### Stainless steel litter pans\*

611536 For 3 sq ft cage  
610149 For 4 sq ft cage  
611537 For 5 sq ft cage

\* 1/2 Floor and Litter Pan must be used together and replaces standard floor



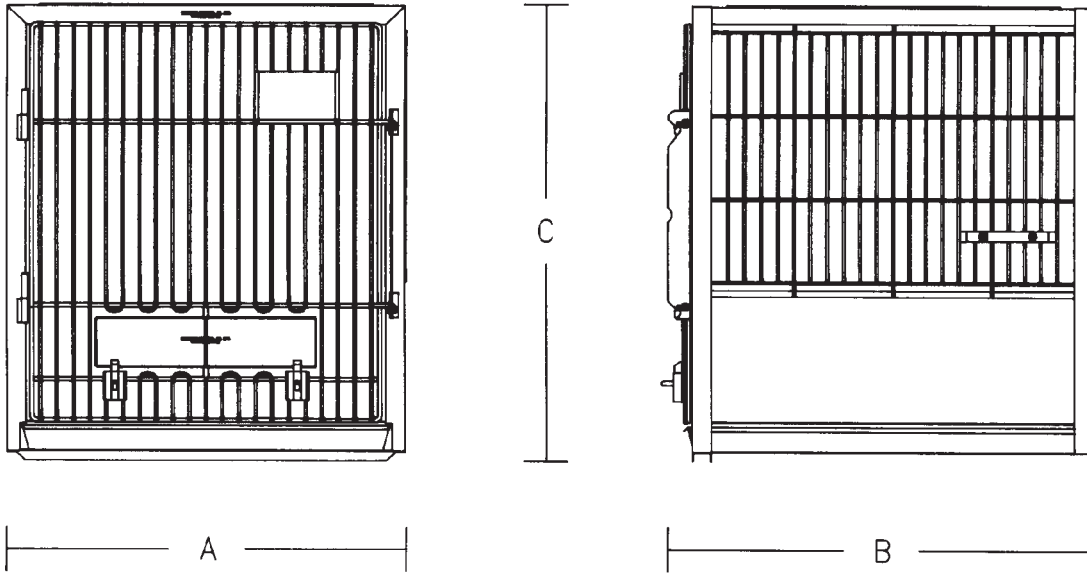

### STAINLESS STEEL DELUXE STACKING CAT CAGES WITH REMOVABLE COLLECTION PANS

PRODUCT NUMBER	ANIMAL SPACE		SHELF*	CAGE DIMENSIONS			FOR ANIMALS WEIGHING**
	FLOOR AREA	ABOVE FLOOR		A	B	C	
101447-00	3.00 sq ft	27.000"	2	18.000"	26.000"	30.875"	≤ 4 kg
101448-00	4.00 sq ft	24.000"	1	24.000"	26.000"	27.422"	> 4 kg
101449-00	5.00 sq ft	27.000"	2	30.000"	26.000"	30.875"	> 4 kg

#### ADDITIONAL EQUIPMENT

Stainless steel rolled edge strip floors	Stainless steel bowls	Stainless steel tapered collection pans
201791 for 3 sq ft cages	200187 .375 qt. bowl	601186 for 3 sq ft cages
201792 for 4 sq ft cages		601187 for 4 sq ft cages
203349 for 5 sq ft cages		603484 for 5 sq ft cages

\* Specify type of shelf: stainless steel, oak, or polyethylene.

\*\* See page on animal space recommendations.

# STAINLESS STEEL ANGULAR RACKS FOR DELUXE STACKING CAT CAGES

Series A Racks support Series A Cages



**101448-00** Cages (shown)  
**104526-00** Rack (shown)

This Stainless Steel Angular Rack for Deluxe Stacking Cat Cages is fabricated from die-formed angles and channels. The angle uprights' lower ends have caster couplers. The angle uprights, angle crossmembers, and channel crossmembers are butted at all joints and welded for strength.

Each cage is held securely in place by a pair of full-length cage runners that are welded to the rack crossmembers.


The rack has a glass beaded to a satin finish, removing any discoloration.

The rack is equipped with chrome plated swivel casters with 5" polyolefin wheels, two with wheel brakes. Swivel and axle bearings have seals and zerk grease fittings.

## OPTIONS

Stainless Steel watering system  
CPVC watering system  
Stainless Steel casters



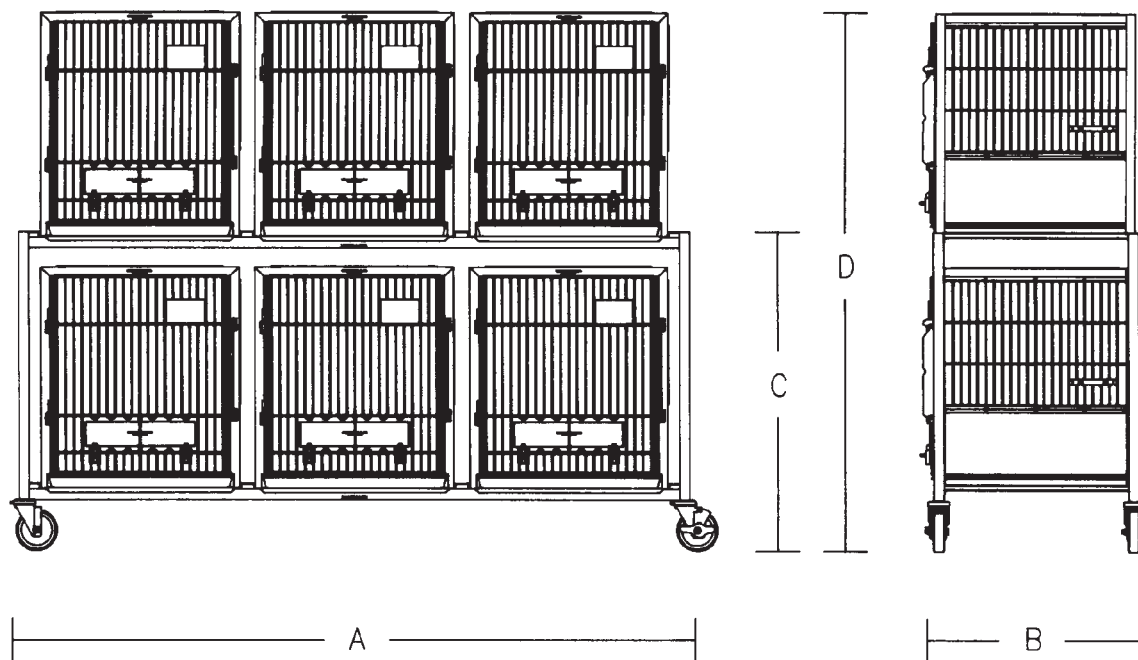

### STAINLESS STEEL ANGULAR RACKS FOR DELUXE STACKING CAT CAGES

PRODUCT NUMBER	RACK CAPACITY			RACK DIMENSIONS			
	SINGLE WIDTH CAGES	CAGES WIDE	TIERS HIGH	A	WITHOUT CAGE B	C	WITH CAGE D
<b>For 3 square foot cages</b>							
104520-00	4	2	2	44.281"	26.625"	40.906"	70.625"
104521-00	6	2	3	44.281"	26.625"	73.656"	103.438"
104522-00	6	3	2	64.188"	26.625"	40.906"	70.625"
104523-00	9	3	3	64.188"	26.625"	73.656"	103.438"
<b>For 4 square foot cages</b>							
104526-00	4	2	2	56.281"	26.625"	37.406"	63.625"
104527-00	6	2	3	56.281"	26.625"	66.656"	92.875"
104528-00	6	3	2	82.188"	26.625"	37.406"	63.625"
104529-00	9	3	3	82.188"	26.625"	66.656"	92.875"
<b>For 5 square foot cages</b>							
104532-00	4	2	2	68.281"	26.625"	40.906"	70.625"
104533-00	6	2	3	68.281"	26.625"	73.656"	103.438"
104534-00	6	3	2	100.188"	26.625"	40.906"	70.625"
104535-00	9	3	3	100.188"	26.625"	73.656"	103.438"

Racks with other cage capacities available

# STAINLESS STEEL DELUXE STACKING METABOLISM CAT CAGES WITH REMOVABLE METABOLISM PANS

Series B Cages fit Series B Racks



**104660-00** Cage (shown)

This Stainless Steel Deluxe Stacking Metabolism Cat Cage with Removable Metabolism Pan features a one-piece wrapper, consisting of one solid side, solid back, and solid lower portion of the opposite side.

The top and upper portion of one side of the cage have electropolished grills. The wrapper and grills are welded to the cage's welded angle frame.

All welds are polished to a satin finish, removing any discoloration.

Each cage is equipped with either one or two shelves, depending on the cage's height. Shelves are available in three types of materials: stainless steel, oak, or polyethylene.

The removable electropolished door has a rod-frame with vertical and horizontal rods, a 5" wide by 3" high cardholder, a rotating lazy-susan feeder with two removable 12 ounce bowls, and a two-point self-locking latch.

The removable electropolished floor with rolled-edge strip animal resting surface, has a frame design that holds the floor off of the wrapper's deflector. This allows waste to drain unobstructed into the removable metabolism pan.

The metabolism pan is supported by the cage's pan runners below the deflector's drip edge. This unique drip edge directs wastes into the metabolism pan. The cage's lower front crossmember is designed to hold the cage securely on the rack.


---

## OPTIONS

---

Stainless steel #2 wire mesh floor\*  
For 3 sq ft cage  
For 4 sq ft cage  
For 5 sq ft cage

\*Mesh floors replace standard floor

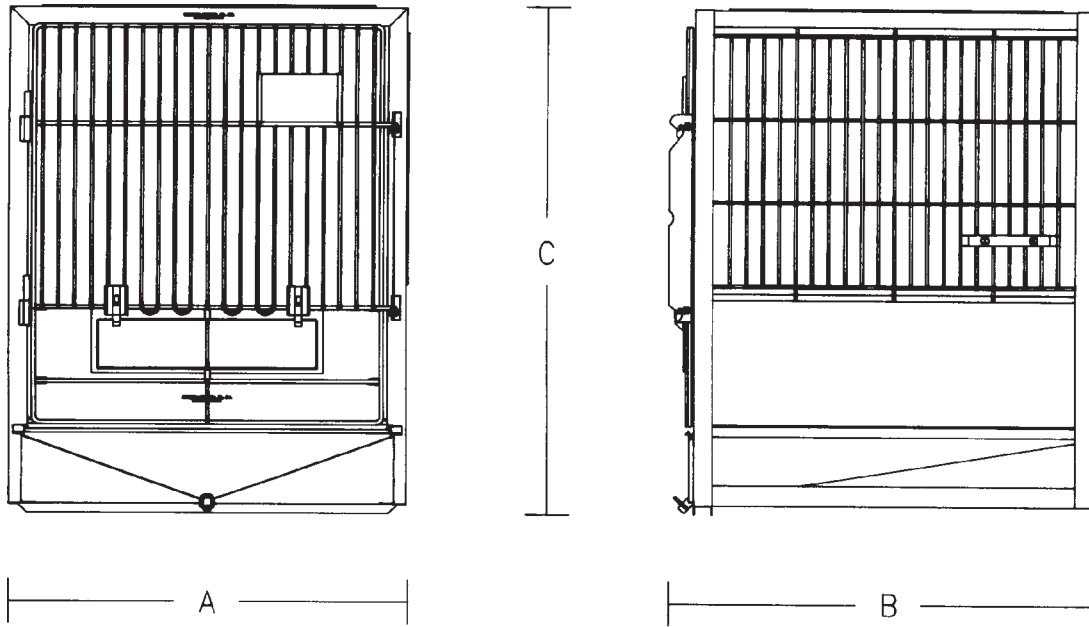

### STAINLESS STEEL DELUXE STACKING CAT CAGES WITH REMOVABLE METABOLISM PANS

PRODUCT NUMBER	ANIMAL SPACE		SHELF*	CAGE DIMENSIONS			FOR ANIMALS WEIGHING**
	FLOOR AREA	ABOVE FLOOR		A	B	C	
104660-00	3.00 sq ft	27.000"	2	18.000"	26.000"	33.437"	< 4 kg
104661-00	4.00 sq ft	24.000"	1	24.000"	26.000"	30.437"	> 4 kg
104662-00	5.00 sq ft	27.000"	2	30.000"	26.000"	33.437"	> 4 kg

#### ADDITIONAL EQUIPMENT


Stainless steel rolled edge strip floors	Stainless steel bowls	Stainless steel metabolism pans
201791 for 3 sq ft cages	200187 .375 qt. bowl	207146 for 3 sq ft cages
201792 for 4 sq ft cages		206774 for 4 sq ft cages
204493 for 5 sq ft cages		207147 for 5 sq ft cages

\* Specify type of shelf: stainless steel , oak, or polypropylene.

\*\* See page on animal recommendations.

# STAINLESS STEEL ANGULAR RACKS FOR DELUXE STACKING METABOLISM CAT CAGES WITH REMOVABLE METABOLISM PANS

Series B Racks support Series B Cages

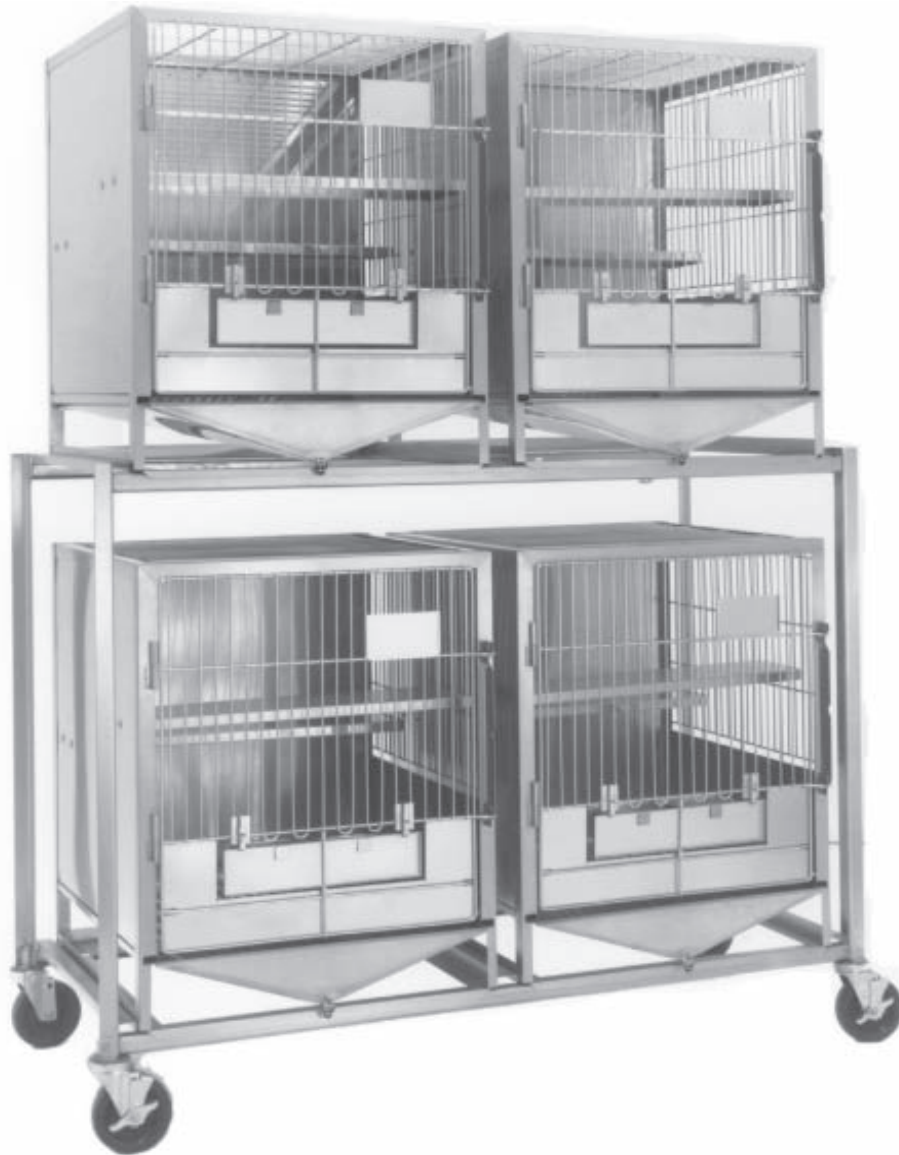

**104661-00** Cages (shown)

**103987-00** Rack (shown)

This Stainless Steel Angular Rack for Deluxe Stacking Metabolism Cat Cages with Removable Metabolism Pans is fabricated from die-formed angles and channels. The angle uprights' lower ends have caster couplers. The angle uprights, angle crossmembers, and channel crossmembers are butted at all joints and welded for strength.

Each cage is held securely in place by a pair of full-length cage runners that are welded to the rack crossmembers.

The rack has a glass beaded satin finish, removing any discoloration.


The rack is equipped with chrome plated swivel casters with 5" polyolefin wheels, two with wheel brakes. Swivel and axle bearings have seals and zerk grease fittings.

---

## OPTIONS

---

Stainless Steel watering system  
CPVC watering system  
Stainless Steel casters

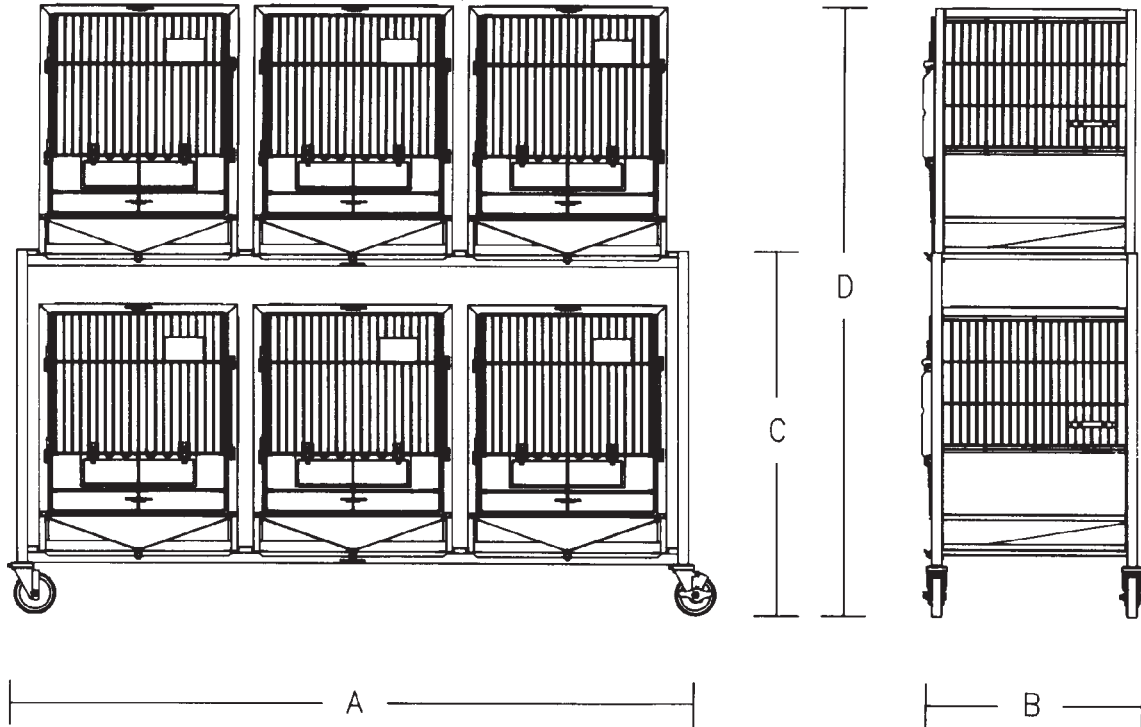

**STAINLESS STEEL ANGULAR RACKS FOR DELUXE STACKING METABOLISM CAT CAGES  
WITH REMOVABLE METABOLISM PANS**

PRODUCT NUMBER	RACK CAPACITY			A	RACK DIMENSIONS		
	SINGLE WIDTH CAGES	CAGES WIDE	TIERS HIGH		WITHOUT CAGE B	C	WITH CAGE D
<b>For 3 square foot cages</b>							
103981-00	4	2	2	44.281"	26.625"	44.375"	73.562"
103982-00	6	2	3	44.281"	26.625"	80.844"	110.031"
103983-00	6	3	2	64.188"	26.625"	44.375"	73.562"
103984-00	9	3	3	64.188"	26.625"	80.844"	110.031"
<b>For 4 square foot cages</b>							
103987-00	4	2	2	56.281"	26.625"	44.375"	73.562"
103988-00	6	2	3	56.281"	26.625"	80.844"	110.031"
103989-00	6	3	2	82.188"	26.625"	44.375"	73.562"
103990-00	9	3	3	82.188"	26.625"	80.844"	110.031"
<b>For 5 square foot cages</b>							
103993-00	4	2	2	68.281"	26.625"	44.375"	73.562"
103994-00	6	2	3	68.281"	26.625"	80.844"	110.031"
103995-00	6	3	2	100.188"	26.625"	44.375"	73.562"
103996-00	9	3	3	100.188"	26.625"	80.844"	110.031"

Racks with other cage capacities available