



STAINLESS STEEL SOCIAL ENVIRONMENT PRIMATE CAGES

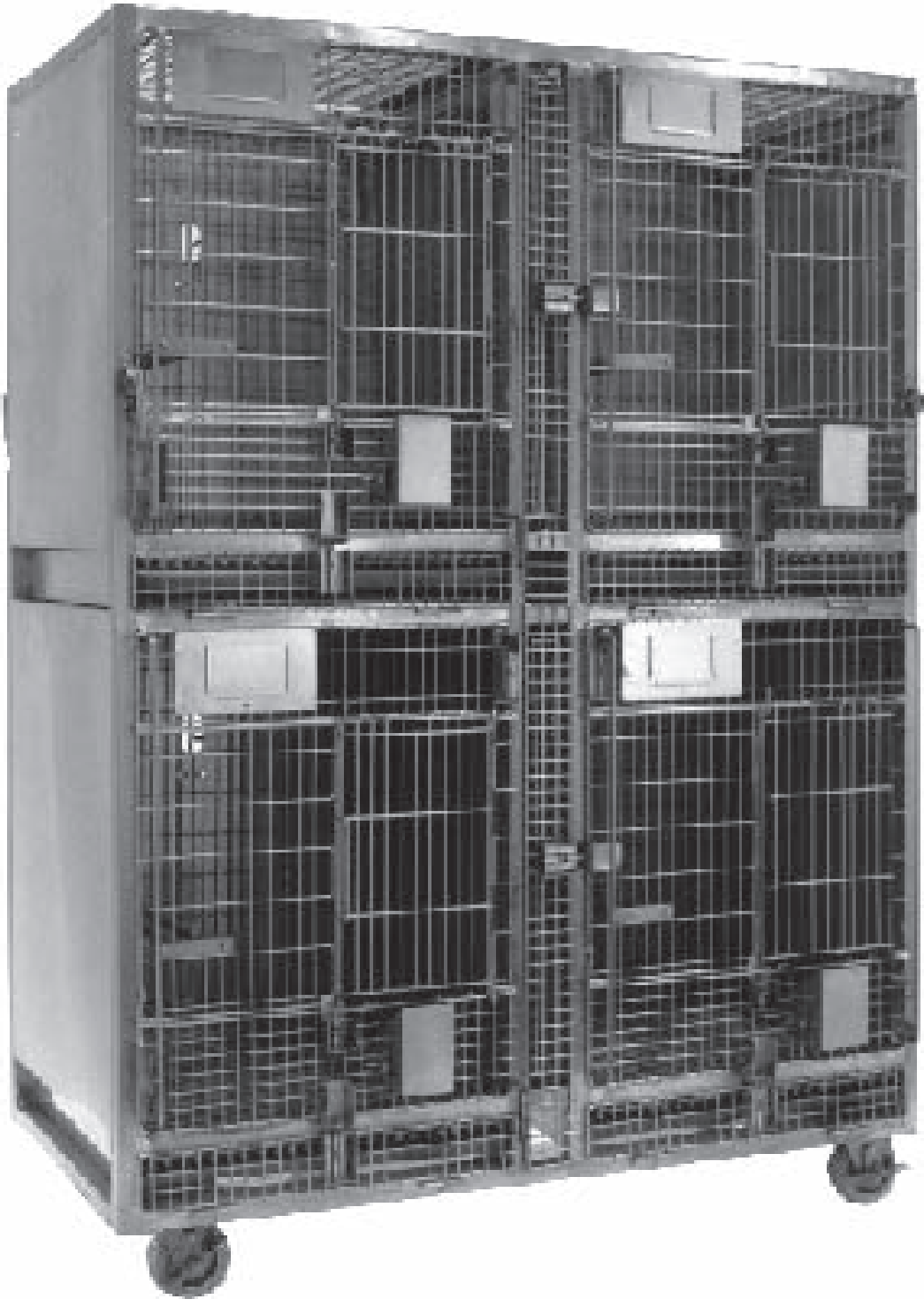

101535-00 Social Environment Primate Cage (Front view shown)

**STAINLESS STEEL SOCIAL ENVIRONMENT PRIMATE CAGES with
OPENING for REMOVABLE SOCIAL ENRICHMENT PANELS and
FLUSH-THROUGH PANS**

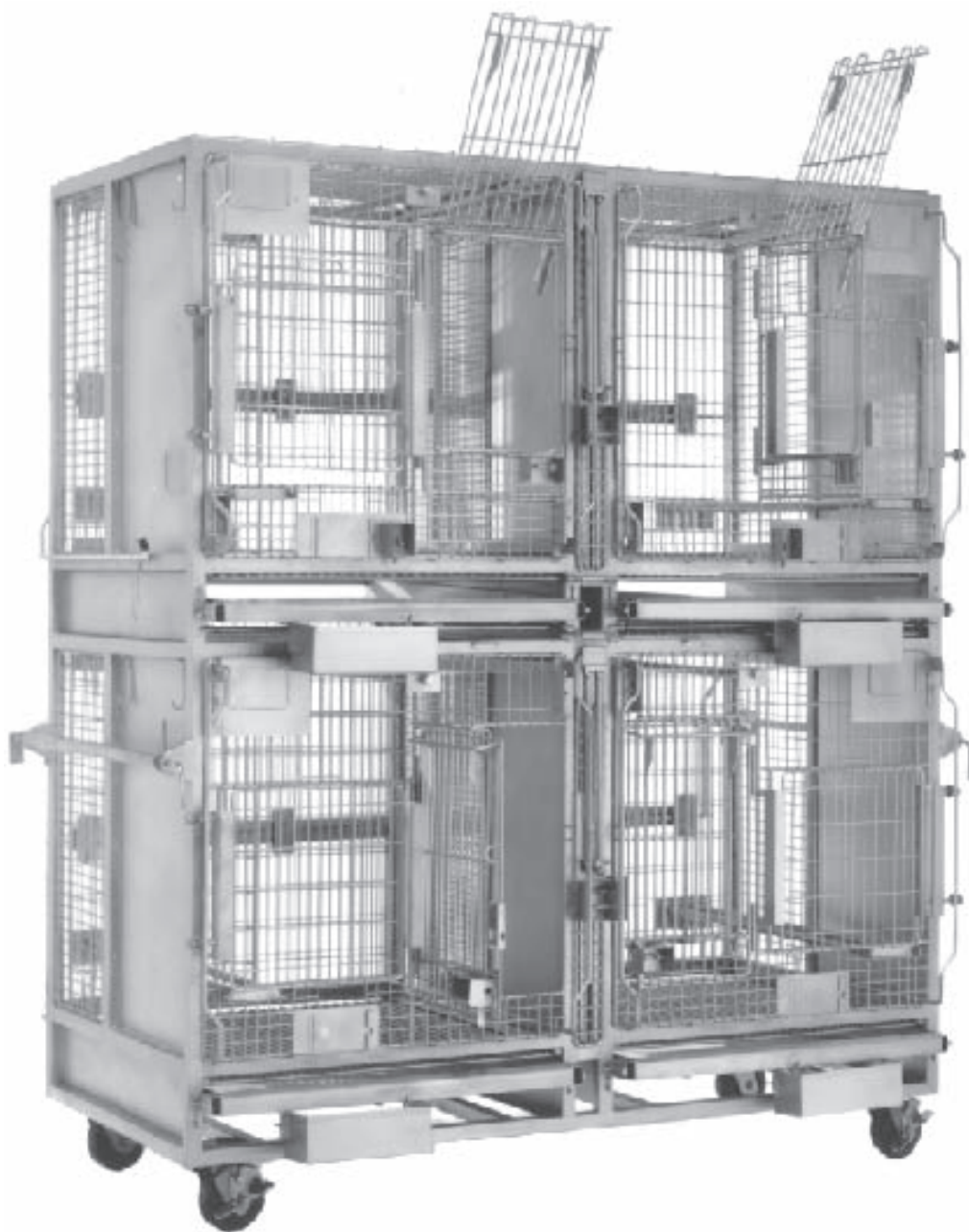

109578 Social Environment Primate Cage (Front view shown)



Opening for Removable Social
Enrichment Panel

**STAINLESS STEEL SOCIAL ENVIRONMENT PRIMATE CAGES with
OPENING for REMOVABLE SOCIAL ENRICHMENT PANELS and
FLUSH-THROUGH PANS**



109578 Social Environment Primate Cage (Back view)

left side has Flush-through Pans removed

right-side has Flush-through Pan in place

STAINLESS STEEL SOCIAL ENVIRONMENT PRIMATE CAGES

This Stainless Steel Social Environment Primate Cage can become a 1, 2, 3, or 4 compartment unit by removing or inserting vertical and horizontal components.

The front of each compartment has fixed rod grill panels fitted with a Guillotine Door, a Sliding Door, or a Hinged Guillotine Door and a 5"W 3"H cardholder, tamper-proof door locks and a removable feeder.

Each compartment's inside features a removable perch and a ratchet-type squeeze back.

Removable floors for the Upper Compartments, removable tops for the Lower Compartments, and removable collection pans for all four compartments are standard. The pans on the upper tier can be removed without the floors/tops being removed.

The removable vertical divider is really 3 dividers: These (3) removable vertical dividers on each level separate the left half from the right. One removable vertical panel has a grill featuring 1" x 1" square mesh, the other has 0.5" x 0.5" square mesh, and the third is fabricated of solid Stainless Steel. Clear acrylic or Stainless Steel with an acrylic "window" panels are available as options. When the vertical dividers are removed, the animals' floor area is more than doubled. A divider with a sliding door is also available, as an option.

The frame is fabricated from 1.25" square tubing. The tubing is butted at all joints and welded for strength. The top, sides, and back can be fitted with wire or solid, removable or fixed panels; in any combination. The frame and all rod grill panels have a glass beaded satin finish, removing any discoloration.

The unit is equipped with chrome plated swivel casters with 5" polyolefin wheels, two with wheel brakes. Swivel and axle bearings have seals and zerk grease fittings.

OPTIONS

Stainless Steel watering system

CPVC watering system

Stainless steel casters

Clear Acrylic Visual Vertical Divider Panel

Stainless Steel Vertical Divider Panel with Clear Acrylic Visual Window

Stainless Steel Vertical Divider Frame with Two Swing Bars

Various Removable Social Enrichment Panel