


# STAINLESS STEEL KEYHOLE SUSPENDED PRIMATE CAGES

Series O Cages fit Series O Racks

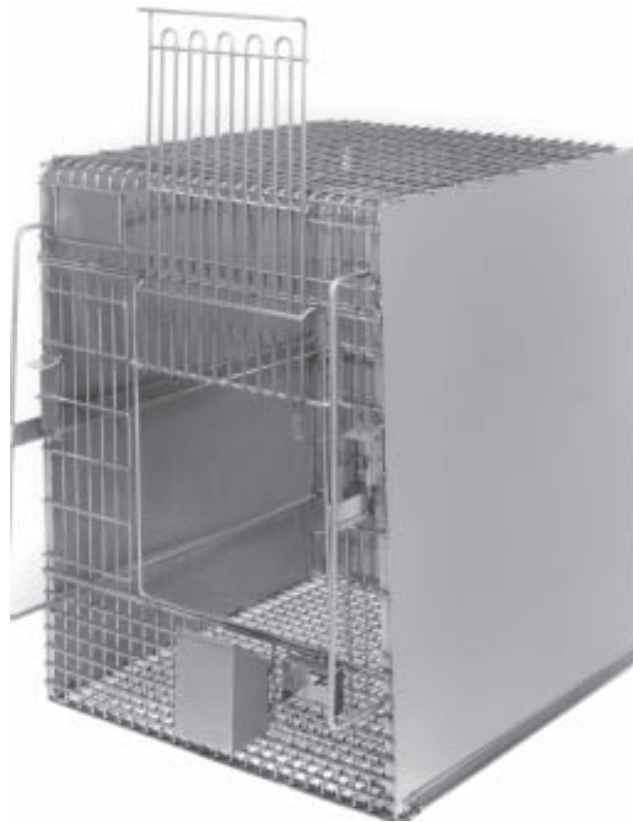

100123-00 Cage (shown)

This Stainless Steel Keyhole Suspended Primate Cage is available with four different sides and back combinations:

Type (1) has a one-piece solid wrapper, consisting of sides and back.

Type (2) has a one-piece solid side and back, with a wire grill on opposite side.

Type (3) has a wire grill on both sides with a solid back.

Type (4) has a one-piece wire wrapper, consisting of sides and back.

The back of each type of cage has full-width reinforcement plates along the top and bottom, which provides additional strength to the back of the cage and an additional area for welding the grill wrapper. The top plate has two keyhole slots, and the bottom has one. The keyhole slots in the reinforcement plate and the keyhole slots in the cage back are fusion welded together around the perimeter. The back of the cage can accommodate a watering system.

All four types of cages have a one-piece grill wrapper, consisting of top, front, and bottom, welded to the sides and back. The front of each cage has a 5" wide by 3" high cardholder. Each cage is equipped with an externally mounted

feeder.

The cage is available with either a guillotine or sliding door; each with a shielded, tamperproof, self-latching catch.

Each cage can be manufactured with or without a ratchet type squeeze back. Also, each cage is designed to support an optional removable metabolism collection pan.

All welds are glass beaded to a satin finish, removing any discoloration.

This cage can be suspended on a wall-mounted system or on various types of racks.

---

## OPTIONS

---



### Stainless Steel Metabolism Collection Pans

|        |                      |
|--------|----------------------|
| 206552 | for 1.6 sq. ft. cage |
| 206553 | for 3.0 sq. ft. cage |
| 206554 | for 4.3 sq. ft. cage |
| 206555 | for 6.0 sq. ft. cage |

Feeder with Back Plate

Sliding Side Door on One (1) or Both Sides



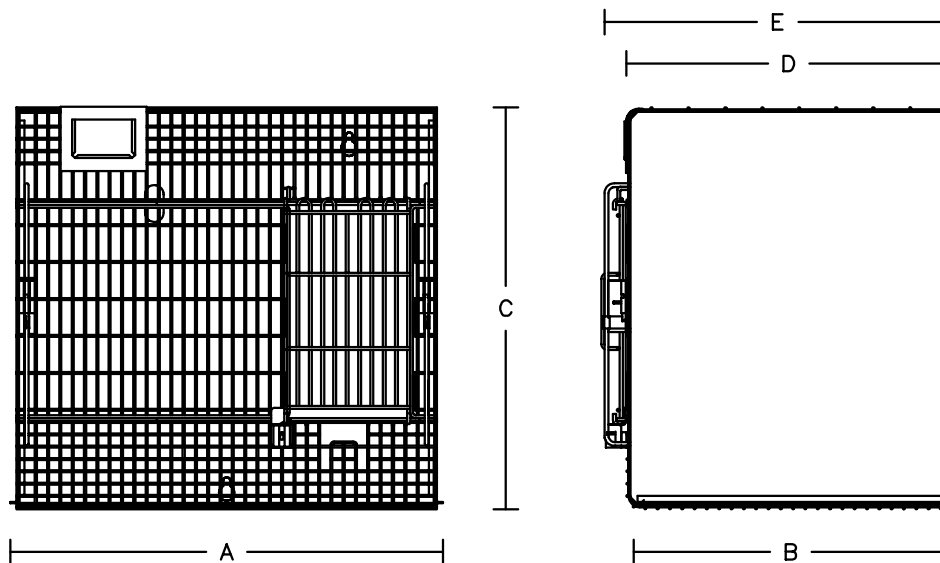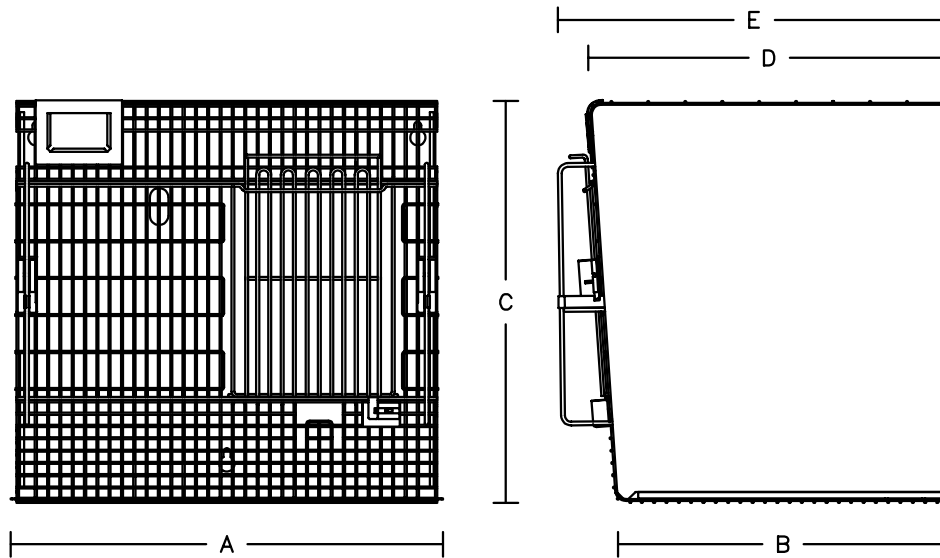



### STAINLESS STEEL KEYHOLE SUSPENDED PRIMATE CAGES\*\*

| PRODUCT NUMBER              | ANIMAL SPACE |             | CAGEDIMENSIONS |         |         |         |         | FOR ANIMALS WEIGHING* |
|-----------------------------|--------------|-------------|----------------|---------|---------|---------|---------|-----------------------|
|                             | FLOOR AREA   | ABOVE FLOOR | A              | B       | C       | D       | E       |                       |
| <b>WITH GUILLOTINE DOOR</b> |              |             |                |         |         |         |         |                       |
| 100005-00                   | 1.60 sq ft   | 22.000"     | 15.048"        | 20.172" | 22.648" | 22.415" | 24.438" | < 1 kg                |
| 100063-00                   | 3.00 sq ft   | 30.000"     | 21.048"        | 27.151" | 30.648" | 29.472" | 31.625" | 1-3 kg                |
| 100123-00                   | 4.30 sq ft   | 30.000"     | 25.098"        | 27.151" | 30.658" | 29.472" | 31.903" | 3-10 kg               |
| 100183-00                   | 6.00 sq ft   | 32.000"     | 35.098"        | 27.151" | 32.642" | 29.472" | 31.903" | 10-15 kg              |
| <b>WITH SLIDING DOOR</b>    |              |             |                |         |         |         |         |                       |
| 100015-00                   | 1.60 sq ft   | 22.000"     | 15.048"        | 20.172" | 22.648" | 20.171" | 24.438" | < 1 kg                |
| 100073-00                   | 3.00 sq ft   | 30.000"     | 21.048"        | 27.151" | 30.648" | 27.151" | 31.625" | 1-3 kg                |
| 100133-00                   | 4.30 sq ft   | 30.000"     | 25.048"        | 27.151" | 30.648" | 27.151" | 31.625" | 3-10 kg               |
| 100193-00                   | 6.00 sq ft   | 32.000"     | 35.048"        | 27.151" | 32.648" | 27.151" | 31.625" | 10-15 kg              |

\* See page on animal space recommendations.

\*\* Specify: type 1, 2, 3, or 4 top, sides, and back combination, and with or without squeeze back

## STAINLESS STEEL KEYHOLE SUSPENDED PRIMATE CAGES MODIFIED WITH OPTIONAL SLIDING SIDE DOOR

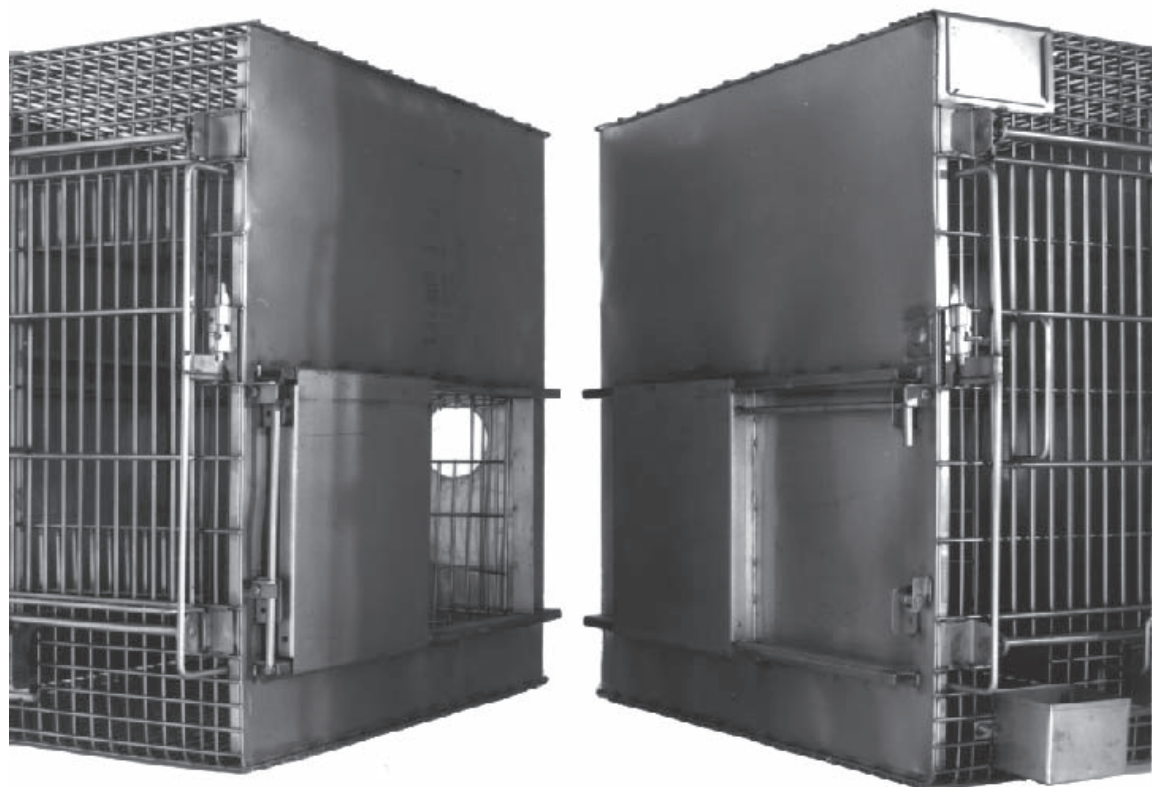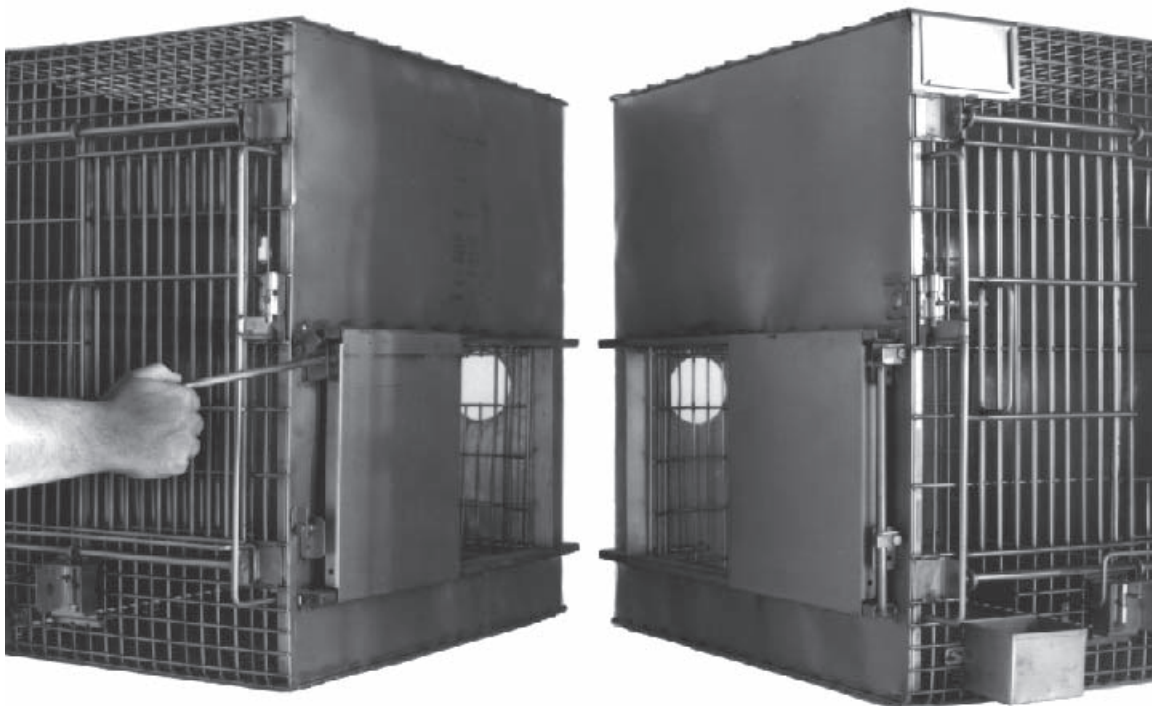


Stainless Steel Suspended Primate Cages can be supplied with an optional Sliding Side Door on one or both sides facilitating socialization.


Your existing cages can also be modified by SSC, Inc. with these doors to facilitate socialization.

[www.subsurg.com](http://www.subsurg.com)

## STAINLESS STEEL VERTICAL SOCIALIZATION TUNNEL FOR PRIMATE CAGES

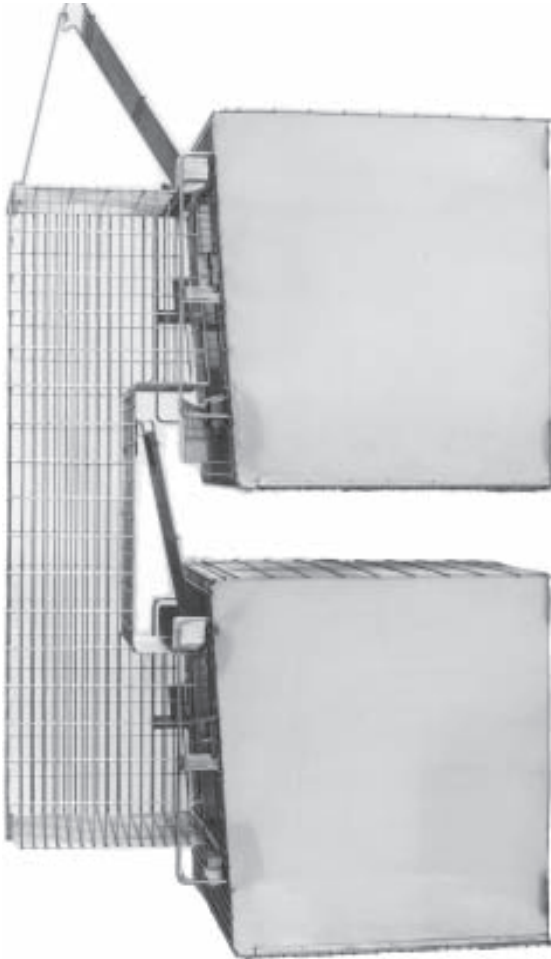

**109285** Tunnel connecting (2) **100123-00** Primate Cages  
(guillotine doors open) shown

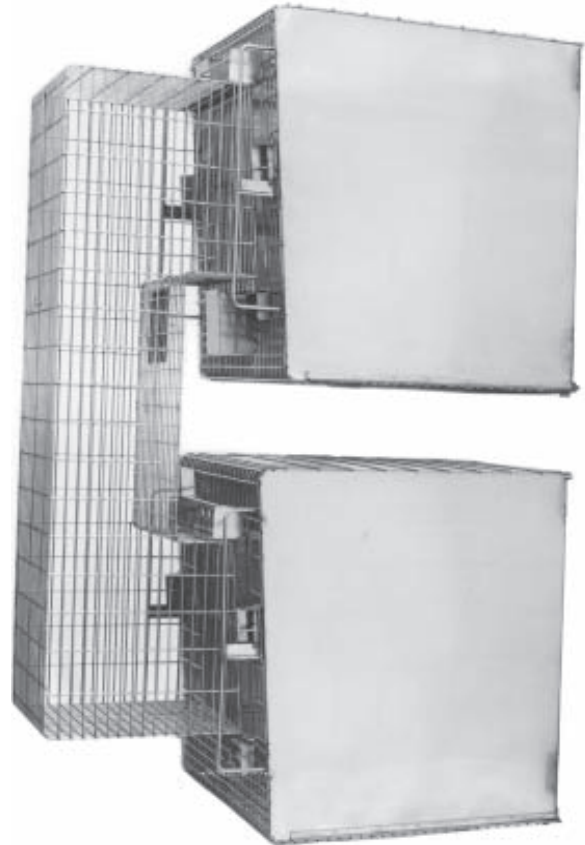

**109285** Tunnel connecting (2) **100123-00** Primate Cages  
(guillotine doors closed) shown

This Vertical Stainless Steel Tunnel for Primate Cages is used to connect two Primate Cages to facilitate socialization.

The Tunnel's design allows the guillotine doors of the cages to operate as normal. Opening the doors allows the animals to interact without touching them.

All welds are polished to a satin finish.

Custom sizes are also available.