

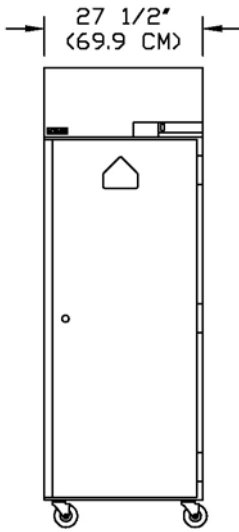


TOP

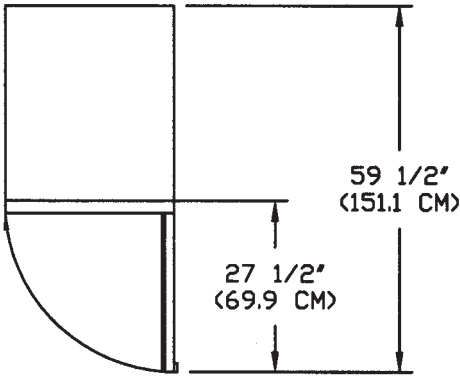

Front

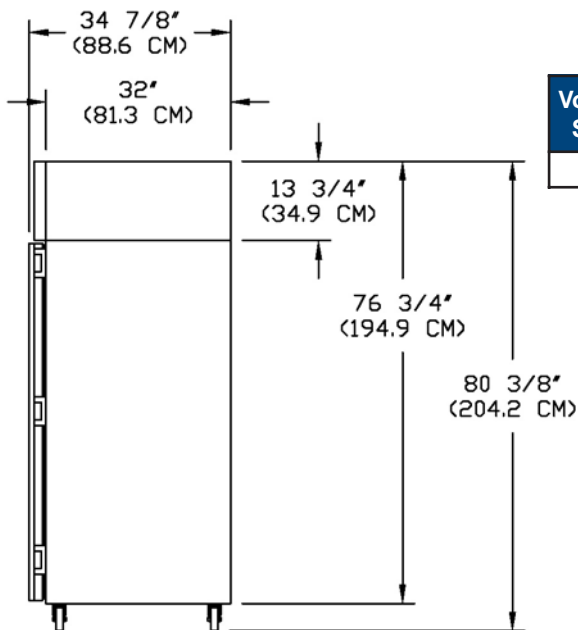

Side

**SPECIFICATIONS**

Model	NSFR241WMW/0M	
Crated Height (in)(cm)	83-3/4"	(212.7 cm)
Crated Width (in)(cm)	37"	(94.0 cm)
Crated Depth (in)(cm)	41"	(104.0 cm)
Crated Weight (lbs)(kg)	325	(148.0 kg)
Interior Height (in)(cm)	59"	(149.9 cm)
Interior Width (in)(cm)	23-1/2"	(59.7 cm)
Interior Depth (in)(cm)	30"	(76.2 cm)
Overall Height* (in)(cm)	80-3/8"	(204.2 cm)
Overall Width (in)(cm)	27-1/2"	(69.9 cm)
Overall Depth (in)(cm)	34-7/8"	(88.6 cm)
Gross Cubage (cu.ft.)(cu.m.)	24	(0.68 cu.m.)
Shelf Area (sq.ft.)(sq.m.)	12.8	(1.19 sq.m.)
Number of Shelves	3	
Number of Casters	4 Standard	
Number of Legs	4 Optional	
Condensing Unit Size	1/4	
Refrigerant	R404A	
Electrical Characteristics	115/60/1	
Cord & Plug Furnished	Yes	
Maximum Fuse Size	15	
NEMA Plug Configuration	5-15P	
Total Amp Draw	6.4	

Specifications subject to change without notice.

Voltage Model Suffix Code	Voltage Description	NEMA Plug	NEMA Receptacle
/0	115V, 1PH, 60HZ	5-15P	5-15R



Revision Date: 7/06  
 © 2006 NOR-LAKE, INC.  
 Printed in the U.S.A.  
 Part Number: 116911



# NOR-LAKE SCIENTIFIC

## FLAMMABLE STORAGE REFRIGERATOR

Designed to meet the demanding requirements of scientific and laboratory applications. Industrial quality cabinet with a heavy duty refrigeration system to maintain your valuable products.


### Standard Product Features:

- Adjustable control range from 2°C to 10°C.  
Factory preset at 4°C
- Internal chamber thermometer
- The exterior front, sides, back and door are white baked on enamel
- The exterior top and bottom are galvanized steel
- Interior liner and interior door panel are smooth aluminum
- Cabinet and doors are foamed-in-place with CFC Free high density polyurethane foam insulation
- Heavy duty spring loaded hinges
- Key lock door
- Full length stainless steel pull handles
- Four casters (2 locking)
- 3 epoxy coated wire shelves
- Manual defrost
- Complete top mount refrigeration system
- Cold wall evaporator
- Air-cooled condensing unit
- Self closing solid door
- Magnetic door gasket
- Warranties: 18 months parts and labor, 5 year compressor(US and Canada) 18 months parts (International)
- UL and C-UL listed; CEC California Energy Commission; Energy Star



### Optional Features:

- Legs in lieu of casters
- Extra shelves
- Stainless steel exterior (with galv. back)
- Stainless steel interior
- Temperature chart recorder
- Chart paper
- Stainless steel drawers (8 per door)

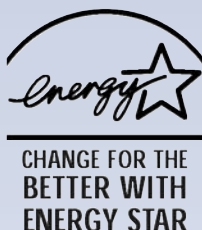

Nor-Lake, Inc.  
Registered to ISO 9001:2000  
File No. A3204

**NOR-LAKE SCIENTIFIC**  
A Division of  
**NOR-LAKE, INCORPORATED**  
727 Second Street  
P.O. Box 248  
Hudson, Wisconsin 54016  
800-477-5253  
715-386-2323  
**800-388-5253 PARTS**  
715-386-4290 FAX  
[www.norlake.com](http://www.norlake.com)