


# NORLAKE MINI-ROOM WALK-INS

SCIENTIFIC®

## TEN DAY SHIPMENT

### Controlled Temperature Chambers-Stability Rooms

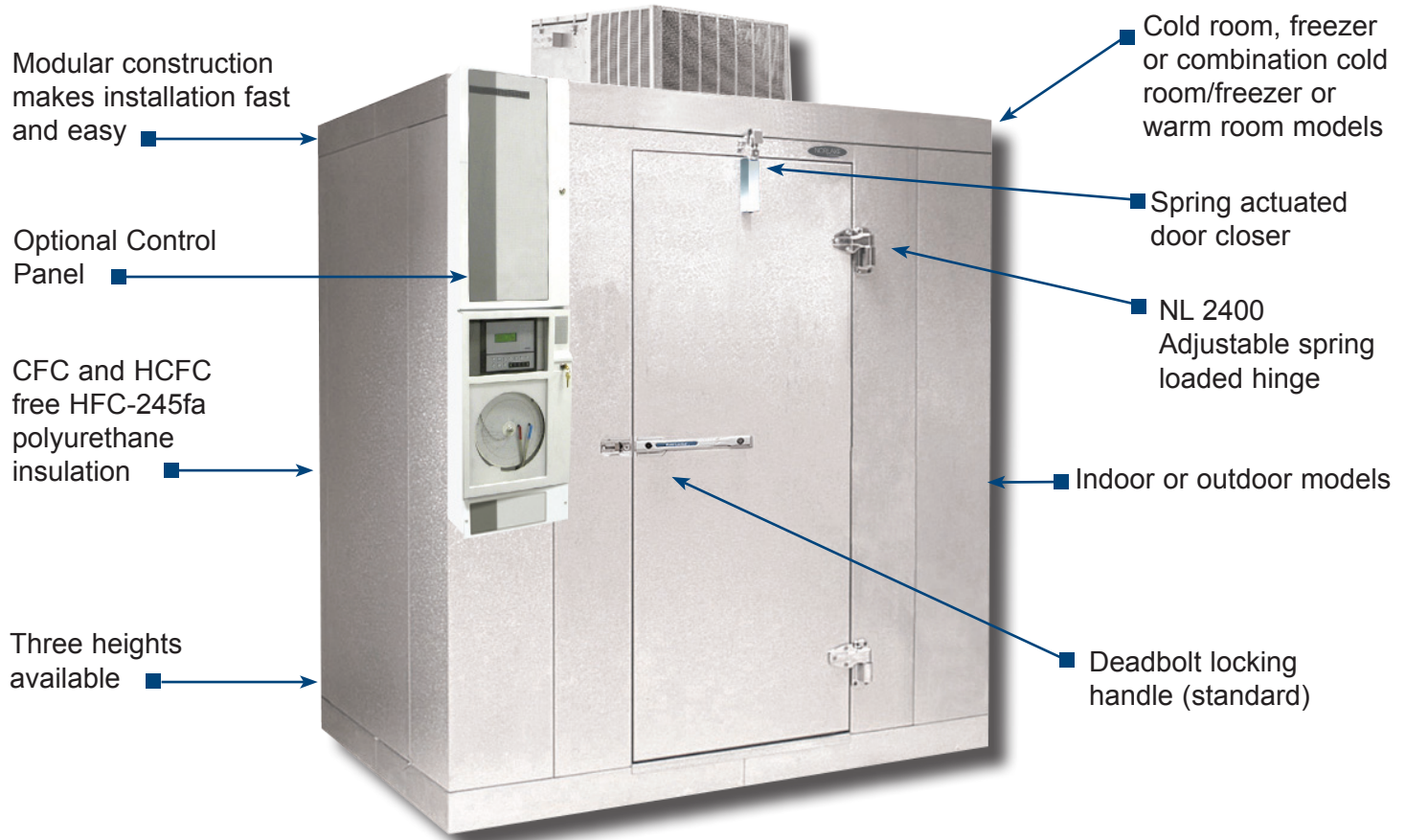

## 10,000 Models Available to Ship in Ten Days!

### Standard Features

- Sizes from 4' x 5' to 8' x 14'
- Three Heights: 6'7" & 7'7" with floor, 7'4" floorless
- Temperatures: +4°C, -20°C, -30°C, +2°C to +55°C
- Exterior finish: 26 ga. corrosion resistant stucco embossed coated steel
- 26" wide self-closing door
- 15 year panel warranty

### Options Available for 10 Day Shipment

- Microprocessor based control panel with accuracies to +0.15°C
- Humidification
- Dehumidification
- Chart recorder
- Viewport door
- Duplex electrical outlets (surface mounted)
- Hi/Lo alarm
- Exterior ramp
- Shelving kits
- Heater kits

### Mini-Room Warm Rooms

- Mini-Room Warm Rooms (+2°C to +55°C)
- Mini-Room Warm Rooms are equipped with a Nor-Lake Scientific Control Panel and a ceiling mount refrigeration system.

**Call 800-477-5253 Today for a Mini-Room Designed Specifically for You.**



## Multiple Applications

- Pharmaceutical Research and Development
- Drug Stability and Shelf Life Testing
- Plasma and Blood Storage
- Clinical Research
- Environment Simulation Testing
- Working Laboratories
- Biological Research
- Electronic Stress Testing


## Mini-Room Walk-in available as...

- Cold room, freezer, combination cold room/freezer or warm room models
- Indoor or outdoor models
- Sizes from 4' x 5' to 8' x 14'
- Three Heights: 6'7" & 7'7" with floor, 7'4" floorless
- Capacity from 77 cubic feet to 1,000 cubic feet
- Wall or ceiling mounted Capsule Pak™ Refrigeration Systems
- Temperatures: +4°C, -20°C, -30°C, +2°C to +55°C



## The Mini-Room Walk-in offers you...

- Installation that is fast and easy
- A full 4" panel foamed-in-place with CFC and HCFC free HFC-245fa polyurethane insulation
- Deadbolt locking handle for security
- Nor-Lake Scientific's exclusive Capsule Pak™ refrigeration system that features a flush coil. 100% usable interior space
- Automatic condensate vaporizer. No plumbing required on indoor units
- Expanded Mini-Room Cold Rooms
  - 6' x 14'
  - 8' x 12'
  - 8' x 14'
  - 8' x 20' combination
- Smaller Mini-Room Walk-ins
  - 3'6" x 6'0" x 6'7" and
  - 3'6" x 7'0" x 6'7"
- Similar floor space as 3 door laboratory reach-in with 20%-30% more capacity

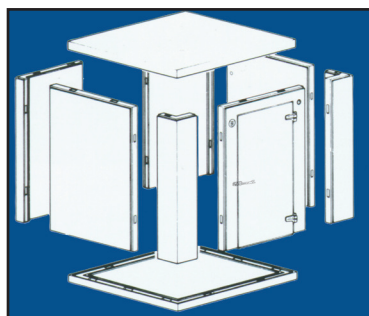

## Modular Assembly

## Capsule Pak™ Refrigeration Systems

- Capsule Pak™ Refrigeration Systems consist of a unitized system which is factory assembled, wired, charged, tested and fully equipped for insertion into a factory prepared walk-in opening. The Capsule Pak™ System has a flush coil which keeps all components outside, allowing more storage inside.
- Wall or ceiling mounted systems
- Indoor or outdoor models
- Temperatures: +4°C, -20°C, -30°C, +2°C to +55°C
- Optional Remote Capsule Pak™ Systems
- UL & C-UL listed

## Shelving Kits

**Specialty designed Shelving Kits are available for all Mini-Room walk-ins.**

- Nor-Lake Scientific Chrome Kote™ Shelving Systems
- Chrome Kote™ shelving available as standard when included
- Quick, easy assembly
- Rigid construction
- Stronger shelves
- Flexible shelf adjustment- 1" increments
- Easy to use
- Nor-Lake shelving utilizes a support clip to eliminate center posts
- The Nor-Lake Klip™ provides maximum space utilization for storage leaving the door opening clear of end standards
- Special epoxy coated finish Green Kote™ shelving provides extra protection from acidic products.

**Contact your sales representative for additional details.**

**NORLAKE** SCIENTIFIC

727 Second Street  
 PO Box 248  
 Hudson, Wisconsin 54016  
 800-477-5253  
 715-386-2323  
 715-386-4290 FAX  
 www.norlakescientific.com  
 © 2014 NOR-LAKE, INC.  
 Revision Date: 06/14  
 Part Number: 089730

