



# Captair Pyramid

Portable glove bag

Safe and convenient transport, storage, and testing of hazardous samples



## Ready-to-use:

- Transport and open suspicious packages
- Sample and user protection
- Protection from potential leaks, or open vials
- Versatility to be used as test chamber

## Portable glove bag

Designed to offer protection from hazardous, known, or potentially contaminated samples. Specimens that are obtained can be placed in the Pyramid, sealed and transported to the lab for testing. If in-field testing is required the Pyramid can be used as a glove box and injected with an inert gas for sample stability.

## Easy to set up and use in the field or lab

It is a multi-function disposable glove bag made of high-quality transparent PVC. Assembled in a few seconds, it can be used anywhere. The slanted shape of the enclosure provides an ergonomic working position for the user.

Light, mobile and disposable, it is an ideal protection tool which can suit many requirements in the laboratory or in the field.


## Optimal air-tightness

Each Pyramid is factory-tested to ensure an airtight seal by inflating with compressed air at a pressure of 2.5 mm Hg.

Once tested, it is given a serial number and a certificate to ensure safety.


## Easy transport

Folds into an easy-to-carry size

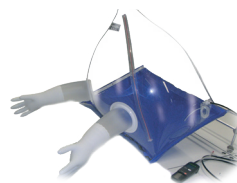

Each unit individually tested and issued a unique serial number for full traceability

# Specifications

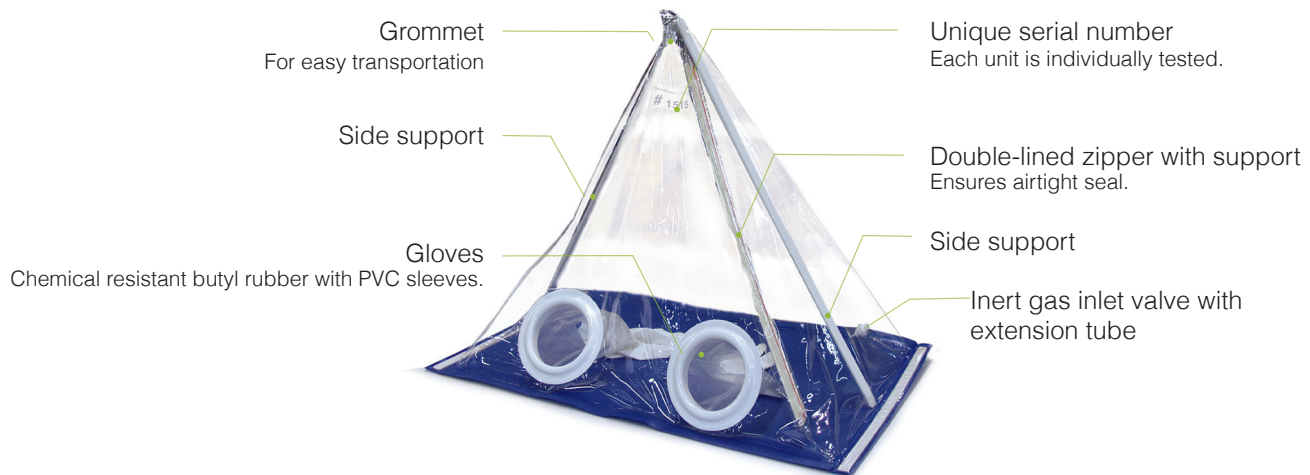

## Dimensions

Width (mm/inches)	Depth (mm/inches)	Height (mm/inches)
860 / 33 <sup>7/8</sup>	560 / 22	725 / 28 <sup>1/2</sup>

## Technical specifications

Enclosure and base	Flexible PVC
Fastener	Double lip zipper (600 mm / 23 <sup>5/8</sup> in)
2 arm openings	Rigid PVC (115 mm / 4 <sup>1/2</sup> in Ø)
Grommet	16 mm / 5/8 in Ø
Operating temperature	-25°C +45°C / -13°F +113°F
Weight	1.55 kg / 3.42 lbs

## Features

Gloves	Chemical resistant butyl rubber with PVC sleeves
Serial number	Unique reference number for better traceability
Valve	Gas insertion valve

**BUY NOW**

For expedited delivery please call 1-800-521-0754 or email [info@labrepco.com](mailto:info@labrepco.com)

**LAB**RepCo

101 Witmer Road, Suite 700, Horsham, PA 19044 Tel: 1-800-521-0754  
[www.labrepco.com](http://www.labrepco.com) [info@labrepco.com](mailto:info@labrepco.com)