



TOP VIEW

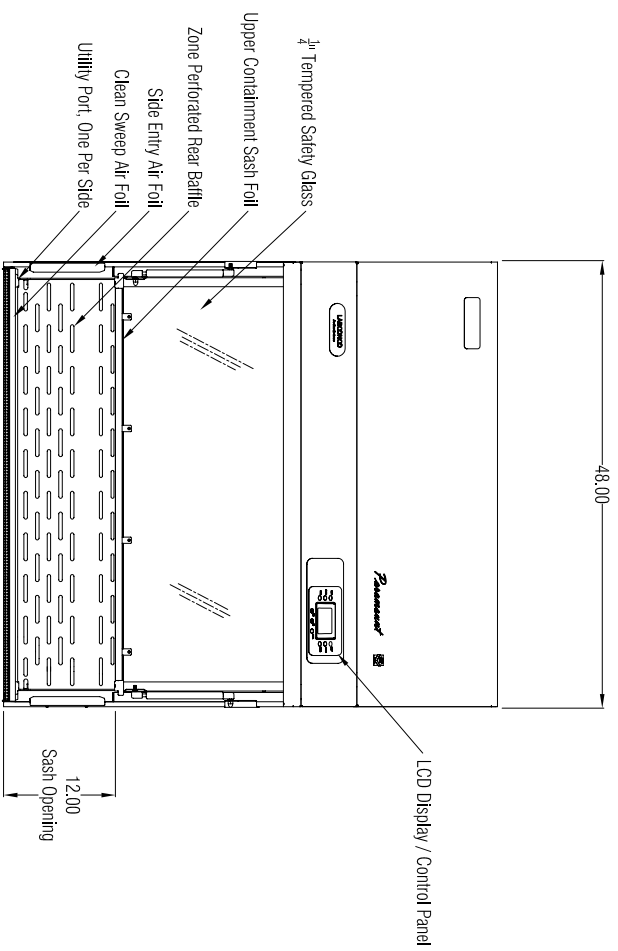

FRONT VIEW

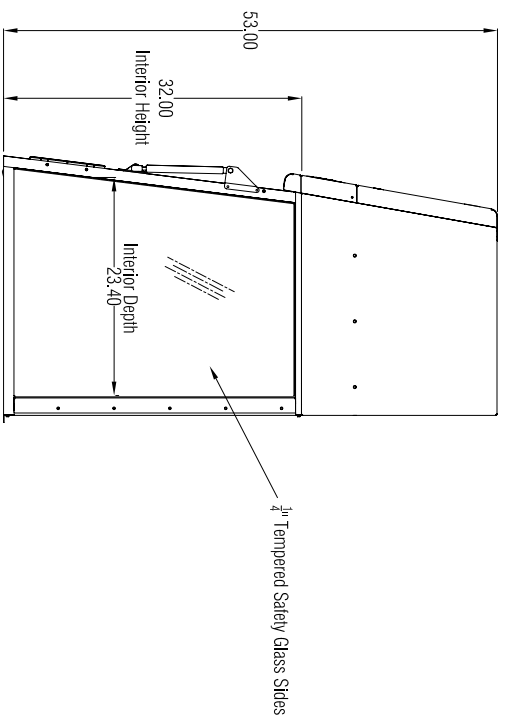

SIDE VIEW

GENERAL DESCRIPTION

4' Paramount Ductless Enclosure, Tall
Model: 6963401

- Exterior Superstructure
 - * Tempered glass
 - * Dry powder epoxy-coated aluminum
- Interior Superstructure
 - * Dry powder epoxy-coated Aluminum and Steel Sides
 - * 1/4" Thick Tempered Safety Glass
- Filters
 - * Accommodates 6 Filters
 - * Gas-assisted filter clamping mechanism

- Electrical Fixtures
 - * Pre-wired Fluorescent Light with Timers

- Control Panel / LCD Display
 - * Face Velocity
 - * Filter Status and Types
 - * Touchpad Control of Enclosure Operation
 - * Set Filter Life and Check Timers
 - * Programmable Alarm Set Points
 - * Event timer

- Wiring
 - * Wired at Factory
 - * 8/2 Power Cord with Plug
- SASH
 - * Gas-Assisted
 - * Pivots to full open for cleaning and material transfer
 - * 1/4" Tempered Safety Glass
 - * 12" Operating Height
 - * Epoxy Coated Aluminum Frame

- Baffles
 - * Perforated, Epoxy Coated, Pivots Downward for Cleaning
 - * Air Foil
 - * Extruded aluminum, Clean Sweep Airfoil
 - Motor
 - * ECM Impeller, Smart-Flow System to Adjust Motor Speed
 - Sensor
 - * Safety First Vapor Sensor
 - * Primary and Secondary Back-Up Sensors
 - Utility Port
 - * (2) 1.5" diameter utility ports with plugs
 - * Assembly and Crating
 - * Pre-Assembled
 - * Palletized
 - Exhaust
 - * Exhaust Shroud, recirculates filtered air back into room

FUME HOODS for:



Labconco Corporation
8811 Prospect Avenue
Kansas City, Mo 64132-2696
(800) 821-5525
www.labconco.com

No.	Date
Date:	4-7-2009
Drawing #	6963401
Drawn By:	Luke Savage
Checked By:	
Approved By:	
Sheet No.	

Order Number: _____ Log Number: _____